



 Sontec Instruments®

LAPAROSCOPY

Modular System

1. Pick a tip
2. Choose a shaft that matches
3. Choose a handle type

**3 items complete
the whole set**


Features of Modular 2 and 3 Piece System

Plastic handle with stepless ratchet

Sontec has fitted this ergonomically-shaped handle with a continuous-action lock that arrests it in any desired position.

A gentle push on the release button is all that is required to unlock it again.

This allows the user to determine the level of pressure to apply to the tissue being clamped, rather than being limited by predefined engaging increments.

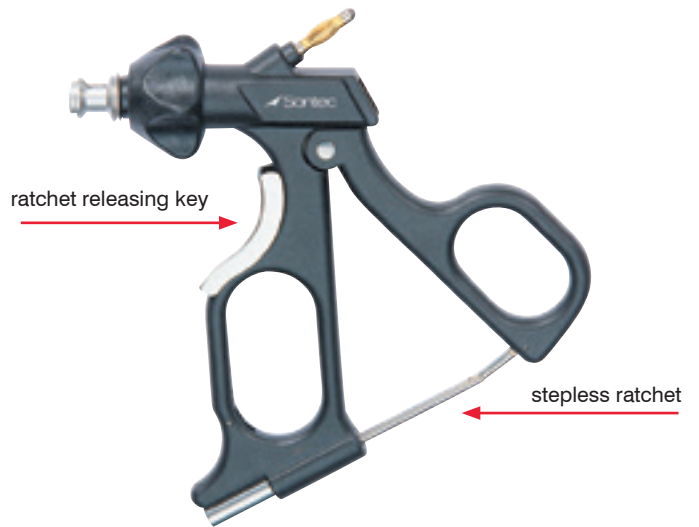

Nut stopper

A nut on the sleeve shaft of Sontec modular instruments prevents bothersome slipping or even the loss of the union nut.




Axial HF connection

The axial HF connection on Sontec modular angled metal handles prevents irritating pivoting motions of the cable. Moreover, the HF connection's location behind the hand aids well-balanced handling of the instrument.

Protective bell fitting for HF connection

A unique protective bell fitting surrounding the HF connection on Sontec modular angled metal handles provides the user with two-fold protection.

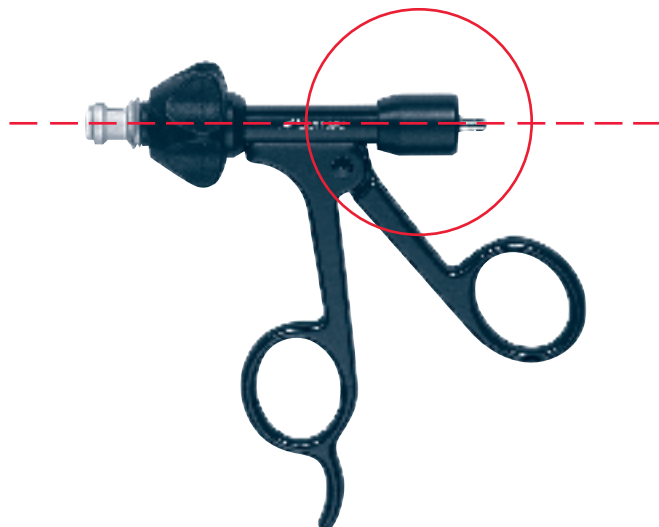


Ergonomic Highly Rigid Plastic Ring Handle

For Ø 3 mm, Ø 5 mm & Ø 10 mm tube shaft instruments

- 2901-090 with stepless ratchet, ratchet activating/release lever, rotation module & HF connection
- 2901-092 with stepless ratchet, ratchet release key, rotation module & HF connection
- 2901-094 with ratchet, ratchet activating/release lever, fold out ratchets, rotation module & HF connection
- 2901-096 with rotation module & HF connection
- 2901-098 with ratchet, ratchet activating/deactivating key, rotation module & HF connection




2901-096

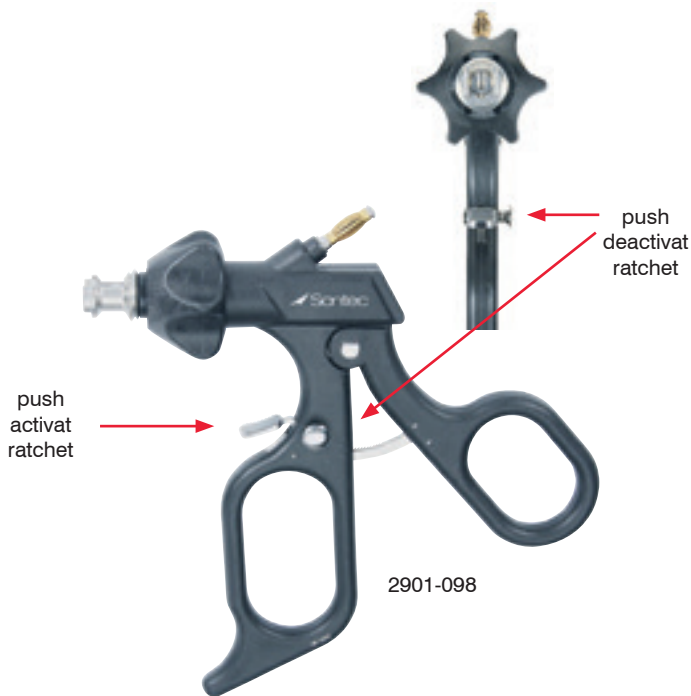

2901-098



2901-092

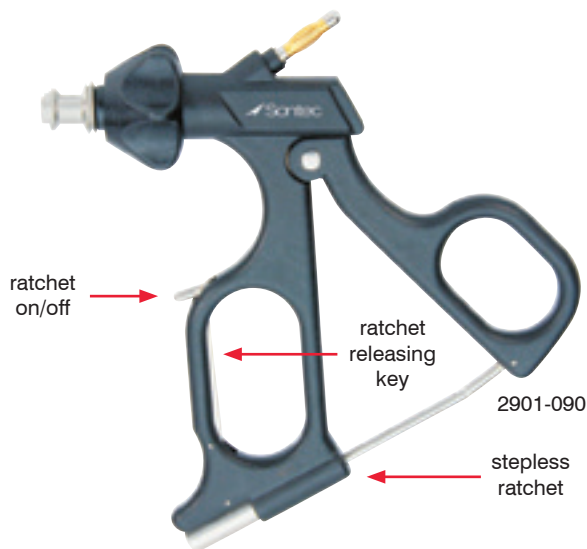

2901-090



2901-094

Angled Insulated Metal Handle

For Ø 3 mm, Ø 5 mm & Ø 10 mm tube shaft instruments

2901-066 with rotating head, ratchet & HF connector



Monopolar Cable

3 meters / 8 mm Bovie Plug /
with 4 mm socket

2695-105



Laparoscopic Bulldog Clamp

Straight

2901-281-45 45 mm jaw
2901-281-60 60 mm jaw



2901-281-45

Curved

2901-281-25 25 mm jaw overall length 2"
2901-281-40 40 mm jaw overall length 2 1/2"



2901-281-60



2901-281-25



2901-281-40

Knot Pusher

With hook / Ø4 mm / 330 mm working length

2901-350



2901-305

Laparoscopic Bulldog Clamp Forceps

Applying clamp / rotatable shaft / articulating jaw / 18" working length

2901-281




2901-281

Axial Ring Handle

For Ø 3 mm, Ø 5 mm & Ø 10 mm tube shaft instruments



- 2901-070 with ratchet
- 2901-072 with ratchet for integral needleholder
- 2901-074 without ratchet
- 2901-078 with branches & ratchet
- 2901-084 with branches & ratchet for integral needleholder angle handle
- 2901-086 with branches & ratchet for integral needleholder angle handle without flushing

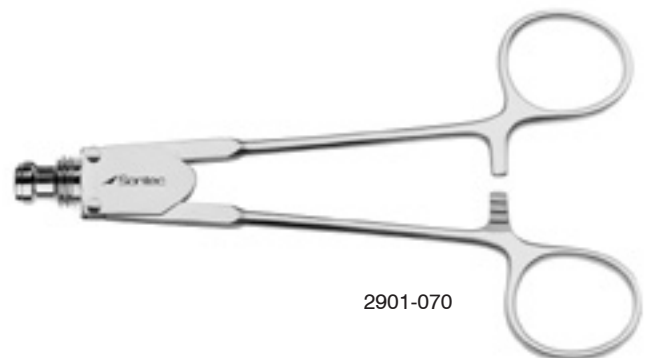

Axial Handle

For Ø 3 mm, Ø 5 mm & Ø 10 mm tube shaft instruments


- 2901-076 with spreading spring & ratchet
- 2901-080 with spreading spring
- 2901-082 with spreading spring & ratchet angled

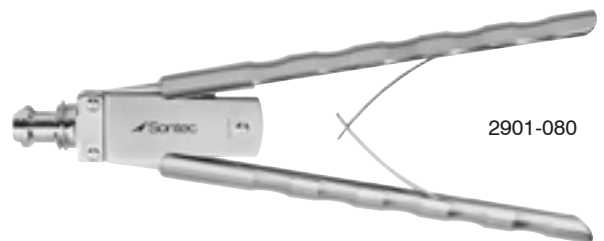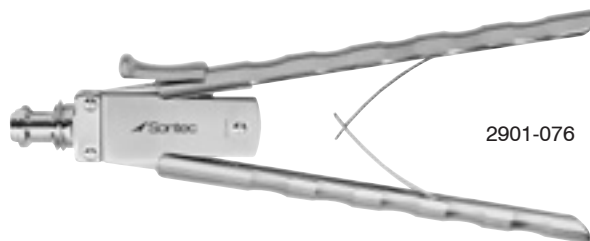

Tube Shaft Insulated

For Ø 5 mm modular 2-piece tube shaft instruments / insulated

2901-101	330 mm shaft length	without LUER lock
2901-101L	330 mm shaft length	with LUER lock
2901-102	360 mm shaft length	without LUER lock
2901-102L	360 mm shaft length	with LUER lock
2901-103	450 mm shaft length	without LUER lock
2901-103L	450 mm shaft length	with LUER lock

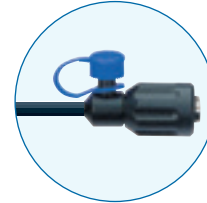

LUER lock



2901-101



2901-102




2901-103

Tube Shaft Insulated

For Ø 10 mm modular 2-piece tube shaft instruments / insulated

2901-104	330 mm shaft length	cylindrical	without LUER lock
2901-104L	330 mm shaft length	cylindrical	with LUER lock
2901-105	330 mm shaft length	reduce to Ø5 mm	without LUER lock
2901-105L	330 mm shaft length	reduce to Ø5 mm	with LUER lock
2901-106	360 mm shaft length	cylindrical	without LUER lock
2901-106L	360 mm shaft length	cylindrical	with LUER lock
2901-107	360 mm shaft length	reduce to Ø5 mm	without LUER lock
2901-107L	360 mm shaft length	reduce to Ø5 mm	with LUER lock
2901-108	450 mm shaft length	cylindrical	without LUER lock
2901-108L	450 mm shaft length	cylindrical	with LUER lock
2901-109	450 mm shaft length	reduce to Ø5 mm	without LUER lock
2901-109L	450 mm shaft length	reduce to Ø5 mm	with LUER lock

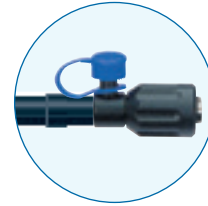

LUER lock



2901-104



2901-105



2901-106



2901-107



2901-108




2901-109

Liver Retractor

Straight triangle / 5-star rotary knob / Ø 5.0 mm


2901-400	60 mm triangle	38 cm length
2901-400L	60 mm triangle	45 cm length
2901-405	80 mm triangle	38 cm length
2901-405L	80 mm triangle	45 cm length
2901-410	120 mm triangle	38 cm length
2901-410L	120 mm triangle	45 cm length

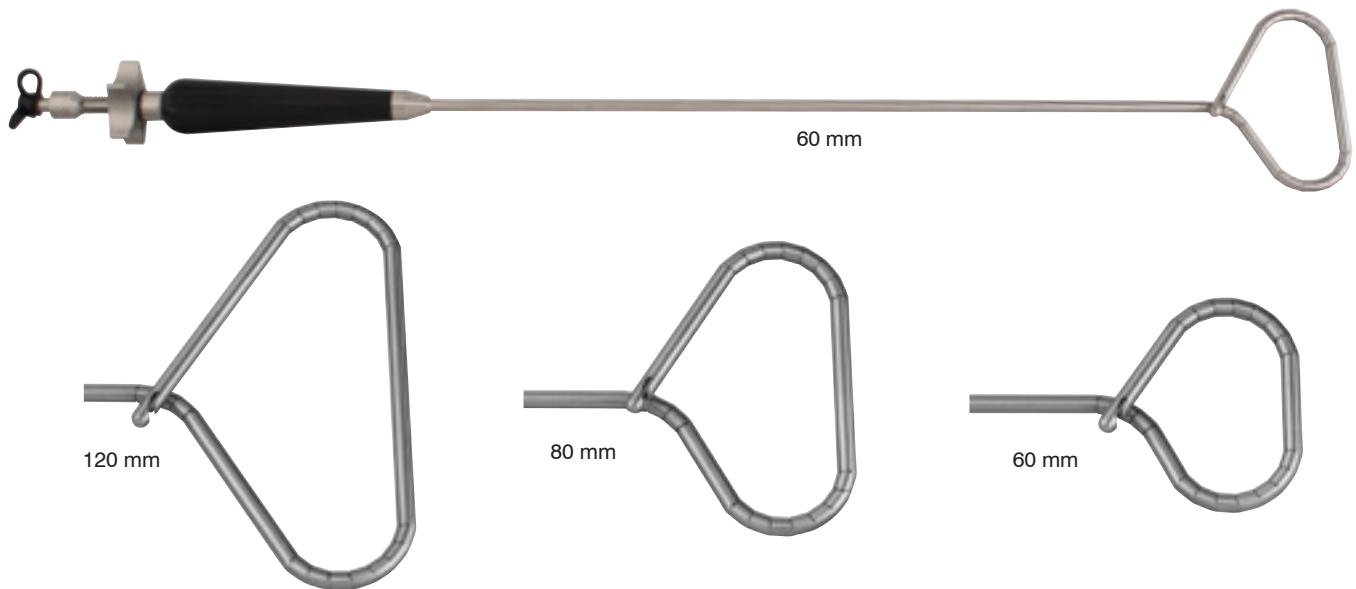

Liver Retractor

Straight triangle / 5-star rotary knob / Ø 5.0 mm

2901-415	60 mm triangle	38 cm length
2901-415L	60 mm triangle	45 cm length
2901-420	80 mm triangle	38 cm length
2901-420L	80 mm triangle	45 cm length

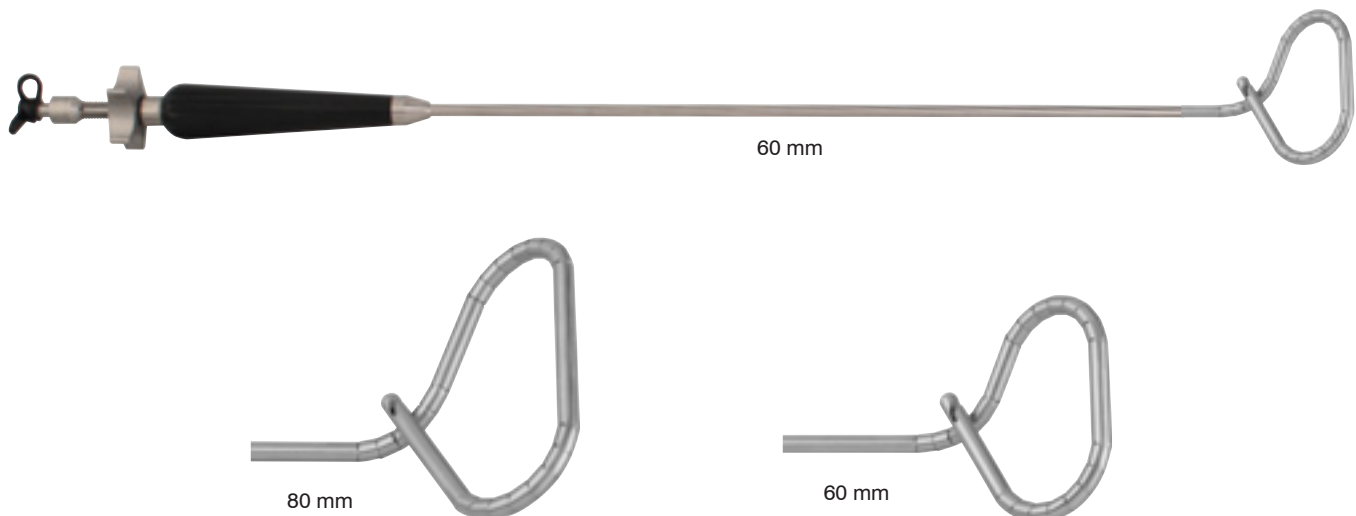

Liver Retractor

Hook / 5-star rotary knob / Ø 5.0 mm

2901-425	90° curved left	38 cm length
2901-425L	90° curved left	45 cm length
2901-430	90° curved right	38 cm length
2901-430L	90° curved right	45 cm length
2901-435	45° curved left	38 cm length
2901-435L	45° curved left	45 cm length
2901-440	45° curved right	38 cm length
2901-440L	45° curved right	45 cm length



90° curved left



90° curved right



45° curved left



45° curved right




45° curved left

Goldfinger Liver Retractor

Finger / Ø 5.0 mm

2901-445-3S	60°	38 cm length	3-star rotary knob
2901-450L-3S	90°	43 cm length	3-star rotary knob
2901-450L	90°t	43 cm length	5-star rotary knob

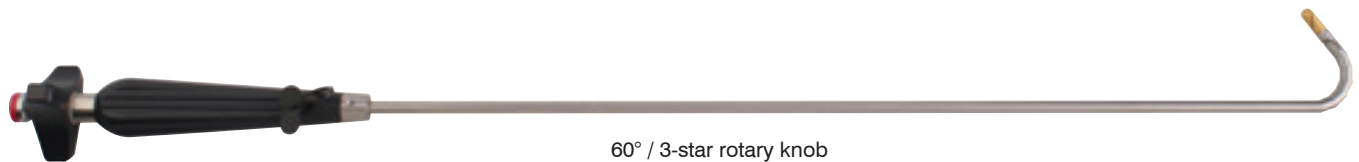

60° / 3-star rotary knob



90° / 5-star rotary knob

Liver Retractor

Round 40 mm / 5-star rotary knob / Ø 5.0 mm

2901-455	straight	38 cm length
2901-455L	straight	45 cm length
2901-460	curved	38 cm length
2901-460L	curved	45 cm length



straight



curved



60 mm

Sontec Fan Organ Retractor

5 blades / 10 mm wide / 350 mm working length

26942-42



26942-42

Esophageal Retractor


1 finger

2901-500 Ø 5 mm 33 cm length

2901-500L Ø 5 mm 45 cm length

2901-501 Ø 10 mm 33 cm length

2901-501L Ø 10 mm 45 cm length

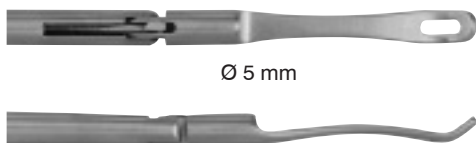

Sontec Fan Organ Retractor

5 blades / 5 mm wide / 330 mm working length / using one hand, blades & neck articulate to achieve desired shape

26942-41



26942-41

Sontec Deluxe Fan Organ Retractor

5 blades / 10 mm wide / 380 mm working length / using one hand, blades & neck articulate to achieve desired shape

26942-43



26942-43

Scissors Insert

For Ø 5 mm tube shaft instruments

Metzenbaum preparation scissors / straight / double action jaw

2901-301S 330 mm shaft length
2901-301 360 mm shaft length
2901-301L 450 mm shaft length



Metzenbaum preparation scissors / curved / double action jaw

2901-302S 330 mm shaft length
2901-302 360 mm shaft length
2901-302L 450 mm shaft length



Power Cut preparation scissors / curved / double action jaw

2901-303S 330 mm shaft length
2901-303 360 mm shaft length
2901-303L 450 mm shaft length



Metzenbaum soft preparation scissors / curved / double action jaw

2901-304S 330 mm shaft length
2901-304 360 mm shaft length
2901-304L 450 mm shaft length



Metzenbaum preparation scissors / curved / double action jaw / serrated / strong version

2901-305S 330 mm shaft length
2901-305 360 mm shaft length
2901-305L 450 mm shaft length



Peritoneal scissors / straight / single action jaw

2901-306S 330 mm shaft length
2901-306 360 mm shaft length
2901-306L 450 mm shaft length



Hook scissors / single action jaw

2901-307S 330 mm shaft length
2901-307 360 mm shaft length
2901-307L 450 mm shaft length



Micro scissors / straight / single action jaw

2901-308S 330 mm shaft length
2901-308 360 mm shaft length
2901-308L 450 mm shaft length



Micro scissors / curved / single action jaw

2901-309S 330 mm shaft length
2901-309 360 mm shaft length
2901-309L 450 mm shaft length



Requires shaft & handle



Grasping Forceps Insert

For Ø 5 mm tube shaft instruments


Universal grasping forceps / double action jaw

2901-211S 330 mm shaft length
2901-211 360 mm shaft length
2901-211L 450 mm shaft length



Maxigrip grasping forceps / double action jaw

2901-212S 330 mm shaft length
2901-212 360 mm shaft length
2901-212L 450 mm shaft length

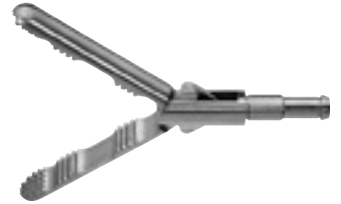

Maryland grasping forceps / double action jaw

2901-213S 330 mm shaft length
2901-213 360 mm shaft length
2901-213L 450 mm shaft length



Kelly grasping forceps / double action jaw

2901-214S 330 mm shaft length
2901-214 360 mm shaft length
2901-214L 450 mm shaft length



Dissecting forceps / fenestrated jaw / double action jaw

2901-215S 330 mm shaft length
2901-215 360 mm shaft length
2901-215L 450 mm shaft length



Dorsey grasping forceps / fenestrated / 20 mm jaw length / double action jaw

2901-216S 330 mm shaft length
2901-216 360 mm shaft length
2901-216L 450 mm shaft length



Dorsey grasping forceps / fenestrated / 20 mm jaw length / single action jaw

2901-217S 330 mm shaft length
2901-217 360 mm shaft length
2901-217L 450 mm shaft length



Dorsey grasping forceps / fenestrated / 40 mm jaw length / double action jaw

2901-218S 330 mm shaft length
2901-218 360 mm shaft length
2901-218L 450 mm shaft length



Duck-Bill grasping forceps / fenestrated / single action jaw

2901-219S 330 mm shaft length
2901-219 360 mm shaft length
2901-219L 450 mm shaft length



Grasping Forceps Insert

For Ø 5 mm tube shaft instruments

Grasping forceps / 2x4 teeth / double action jaw

2901-220S 330 mm shaft length
2901-220 360 mm shaft length
2901-220L 450 mm shaft length



Grasping forceps / 2x4 teeth / single action jaw

2901-221S 330 mm shaft length
2901-221 360 mm shaft length
2901-221L 450 mm shaft length



Grasping forceps / with reservoir / double action jaw / blunt

2901-222S 330 mm shaft length
2901-222 360 mm shaft length
2901-222L 450 mm shaft length



Grasping forceps / with reservoir / double action jaw / sharp

2901-223S 330 mm shaft length
2901-223 360 mm shaft length
2901-223L 450 mm shaft length



Grasping forceps / double action jaw

2901-224S 330 mm shaft length
2901-224 360 mm shaft length
2901-224L 450 mm shaft length



Dolphin grasping forceps / double action jaw

2901-225S 330 mm shaft length
2901-225 360 mm shaft length
2901-225L 450 mm shaft length



Dolphin grasping forceps / long / with reservoir / double action jaw

2901-226S 330 mm shaft length
2901-226 360 mm shaft length
2901-226L 450 mm shaft length



Grasping forceps / smooth / double action jaw

2901-227S 330 mm shaft length
2901-227 360 mm shaft length
2901-227L 450 mm shaft length



Grasping forceps / sharp / crosswise teeth / 14 mm jaw length / double action jaw

2901-228S 330 mm shaft length
2901-228 360 mm shaft length
2901-228L 450 mm shaft length



Requires shaft & handle



Grasping Forceps Insert

For Ø 5 mm tube shaft instruments

Grasping forceps / single action jaw / straight / crosswise teeth / 18 mm jaw length

2901-229S 330 mm shaft length
 2901-229 360 mm shaft length
 2901-229L 450 mm shaft length



Grasping forceps / double action jaw / sharp / crosswise teeth / 18 mm jaw length

2901-230S 330 mm shaft length
 2901-230 360 mm shaft length
 2901-230L 450 mm shaft length



Grasping forceps / distal / with row of teeth / single action jaw

2901-231S 330 mm shaft length
 2901-231 360 mm shaft length
 2901-231L 450 mm shaft length



Grasping forceps / distal / with row of teeth / double action jaw

2901-232S 330 mm shaft length
 2901-232 360 mm shaft length
 2901-232L 450 mm shaft length



Grasping forceps / double action jaw

2901-233S 330 mm shaft length
 2901-233 360 mm shaft length
 2901-233L 450 mm shaft length



Grasping forceps / with fine crosswise teeth / double action jaw

2901-234S 330 mm shaft length
 2901-234 360 mm shaft length
 2901-234L 450 mm shaft length



Grasping forceps / with mouse teeth / double action jaw

2901-235S 330 mm shaft length
 2901-235 360 mm shaft length
 2901-235L 450 mm shaft length



Grasping forceps / with mouse teeth / single action jaw

2901-236S 330 mm shaft length
 2901-236 360 mm shaft length
 2901-236L 450 mm shaft length



Requires shaft & handle



Grasping Forceps Insert

For Ø 5 mm tube shaft instruments

Grasping forceps / with DeBakey teeth / double action jaw

2901-237S 330 mm shaft length
2901-237 360 mm shaft length
2901-237L 450 mm shaft length



Extracting forceps / with rows of retrograde teeth / double action jaw

2901-238S 330 mm shaft length
2901-238 360 mm shaft length
2901-238L 450 mm shaft length



Mixer grasping forceps / 90° angle / double action jaw

2901-239S 330 mm shaft length
2901-239 360 mm shaft length
2901-239L 450 mm shaft length



Mixer grasping forceps / 45° angle / double action jaw

2901-240S 330 mm shaft length
2901-240 360 mm shaft length
2901-240L 450 mm shaft length



Tenaculum grasping forceps / double action jaw

2901-241S 330 mm shaft length
2901-241 360 mm shaft length
2901-241L 450 mm shaft length




Claw grasping & extracting forceps / 2x3 teeth / single action jaw

2901-242S 330 mm shaft length
2901-242 360 mm shaft length
2901-242L 450 mm shaft length



AGRESSIV claw grasping forceps / double action jaw

2901-243S 330 mm shaft length
2901-243 360 mm shaft length
2901-243L 450 mm shaft length

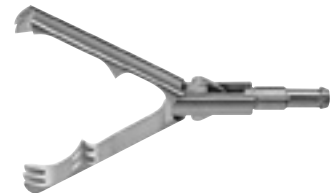

Requires shaft & handle



Grasping Forceps Insert

For Ø 5 mm tube shaft instruments


Atraumatic grasping forceps / with wavy teeth / with reservoir / single action jaw

2901-244S 330 mm shaft length
2901-244 360 mm shaft length
2901-244L 450 mm shaft length



Atraumatic grasping forceps / with wavy teeth / with reservoir / double action jaw

2901-245S 330 mm shaft length
2901-245 360 mm shaft length
2901-245L 450 mm shaft length

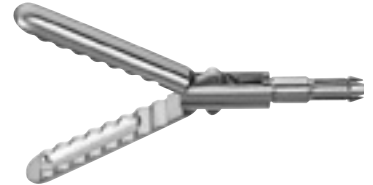

Claw grasping forceps for retraction / with wavy teeth / single action jaw

2901-246S 330 mm shaft length
2901-246 360 mm shaft length
2901-246L 450 mm shaft length




Claw grasping forceps for retraction / with wavy teeth / double action jaw

2901-247S 330 mm shaft length
2901-247 360 mm shaft length
2901-247L 450 mm shaft length



WAVE grasping forceps / fenestrated / double action jaw

2901-248S 330 mm shaft length
2901-248 360 mm shaft length
2901-248L 450 mm shaft length

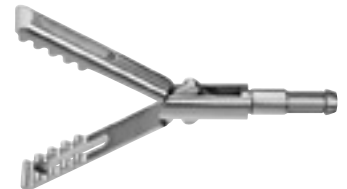

Clamp Insert

For Ø 5 mm tube shaft instruments

DeBakey clamp / 30 mm jaw length / double action jaw

2901-249S 330 mm shaft length
2901-249 360 mm shaft length
2901-249L 450 mm shaft length



Babcock clamp / double action jaw

2901-250S 330 mm shaft length
2901-250 360 mm shaft length
2901-250L 450 mm shaft length



Cup Insert

For Ø 5 mm tube shaft instruments

Tumor grasping forceps / pointed / single action jaw

2901-121S 330 mm shaft length
2901-121 360 mm shaft length
2901-121L 450 mm shaft length



Biopsy cup / double action jaw

2901-123S 330 mm shaft length
2901-123 360 mm shaft length
2901-123L 450 mm shaft length



Biopsy cup / with thorn / double action jaw

2901-125S 330 mm shaft length
2901-125 360 mm shaft length
2901-125L 450 mm shaft length



Biopsy cup long / single action jaw

2901-127S 330 mm shaft length
2901-127 360 mm shaft length
2901-127L 450 mm shaft length



Biopsy cup long / with thorn / single action jaw

2901-129S 330 mm shaft length
2901-129 360 mm shaft length
2901-129L 450 mm shaft length



Alligator grasping forceps / double action jaw

2901-131S 330 mm shaft length
2901-131 360 mm shaft length
2901-131L 450 mm shaft length



Punch Insert

For Ø5 mm tube shaft instruments

Biopsy punch / without thorn / single action jaw

2901-115S 330 mm shaft length
2901-115 360 mm shaft length
2901-115L 450 mm shaft length



Biopsy punch / with thorn / single action jaw

2901-117S 330 mm shaft length
2901-117 360 mm shaft length
2901-117L 450 mm shaft length



Requires shaft & handle



Needle Holder Insert

For Ø 5 mm tube shaft instruments

Needle holder / with TC serrated crosswise teeth / straight / single action jaw

2901-321S 330 mm shaft length
2901-321 360 mm shaft length
2901-321L 450 mm shaft length



Needle holder / with TC serrated crosswise teeth / left curved / single action jaw

2901-322S 330 mm shaft length
2901-322 360 mm shaft length
2901-322L 450 mm shaft length



Needle holder / with TC serrated crosswise teeth / right curved / single action jaw

2901-323S 330 mm shaft length
2901-323 360 mm shaft length
2901-323L 450 mm shaft length



Needle holder positioner / with TC dust / convex/concave shaped jaws with fixation teeth / single action jaw

2901-324S 330 mm shaft length
2901-324 360 mm shaft length
2901-324L 450 mm shaft length



Dissecting Forceps Insert

For Ø5 mm tube shaft instruments

Dissecting forceps / fenestrated / double action jaw

2901-252S 330 mm shaft length
2901-252 360 mm shaft length
2901-252L 450 mm shaft length



Dissecting forceps / fenestrated / slim / double action jaw

2901-254S 330 mm shaft length
2901-254 360 mm shaft length
2901-254L 450 mm shaft length



Requires shaft & handle



Grasping Forceps Insert

For Ø 10 mm tube shaft instruments

Claw grasping & extracting forceps / 2x3 teeth / single action jaw

2901-260S 330 mm shaft length
2901-260 360 mm shaft length
2901-260L 450 mm shaft length



Claw grasping & extracting forceps / 2x3 teeth / double action jaw

2901-261S 330 mm shaft length
2901-261 360 mm shaft length
2901-261L 450 mm shaft length



Maryland grasping forceps / double action jaw

2901-262S 330 mm shaft length
2901-262 360 mm shaft length
2901-262L 450 mm shaft length



Grasping forceps / with DeBakey teeth / fenestrated / double action jaw

2901-264S 330 mm shaft length
2901-264 360 mm shaft length
2901-264L 450 mm shaft length



Babcock grasping forceps / double action jaw

2901-265S 330 mm shaft length
2901-265 360 mm shaft length
2901-265L 450 mm shaft length



Babcock grasping forceps / with crosswise teeth / double action jaw

2901-266S 330 mm shaft length
2901-266 360 mm shaft length
2901-266L 450 mm shaft length



Requires shaft & handle



Grasping Forceps Insert

For Ø 10 mm tube shaft instruments

Mixer grasping forceps / double action jaw

2901-267S 330 mm shaft length
2901-267 360 mm shaft length
2901-267L 450 mm shaft length



Dorsey grasping forceps / double action jaw

2901-268S 330 mm shaft length
2901-268 360 mm shaft length
2901-268L 450 mm shaft length



Tenaculum forceps / double action jaw

2901-275S 330 mm shaft length
2901-275 360 mm shaft length
2901-275L 450 mm shaft length



Cup Insert

For Ø 10 mm tube shaft instruments

Large cup / single action jaw

2901-280S 330 mm shaft length
2901-280 360 mm shaft length
2901-280L 450 mm shaft length



Scissors Insert

For Ø 10 mm tube shaft instruments

Large operating scissors / single action jaw

2901-279S 330 mm shaft length
2901-279 360 mm shaft length
2901-279L 450 mm shaft length



Requires shaft & handle



Grasping Clamp Insert

For Ø 10 mm tube shaft instruments

Kocher clamp / double action jaw

2901-269S	30 mm jaw length / 330 mm shaft length
2901-269	30 mm jaw length / 360 mm shaft length
2901-269L	30 mm jaw length / 450 mm shaft length
2901-270S	40 mm jaw length / 330 mm shaft length
2901-270	40 mm jaw length / 360 mm shaft length
2901-270L	40 mm jaw length / 450 mm shaft length
2901-271S	50 mm jaw length / 330 mm shaft length
2901-271	50 mm jaw length / 360 mm shaft length
2901-271L	50 mm jaw length / 450 mm shaft length



DeBakey clamp / double teeth / double action jaw

2901-272S	30 mm jaw length / 330 mm shaft length
2901-272	30 mm jaw length / 360 mm shaft length
2901-272L	30 mm jaw length / 450 mm shaft length
2901-273S	40 mm jaw length / 330 mm shaft length
2901-273	40 mm jaw length / 360 mm shaft length
2901-273L	40 mm jaw length / 450 mm shaft length
2901-274S	50 mm jaw length / 330 mm shaft length
2901-274	50 mm jaw length / 360 mm shaft length
2901-274L	50 mm jaw length / 450 mm shaft length



Babcock clamp / double action jaw

2901-276S	330 mm shaft length
2901-276	360 mm shaft length
2901-276L	450 mm shaft length




Collin clamp / double action jaw

2901-277S	330 mm shaft length
2901-277	360 mm shaft length
2901-277L	450 mm shaft length



Pennington clamp / double action jaw

2901-278S	330 mm shaft length
2901-278	360 mm shaft length
2901-278L	450 mm shaft length

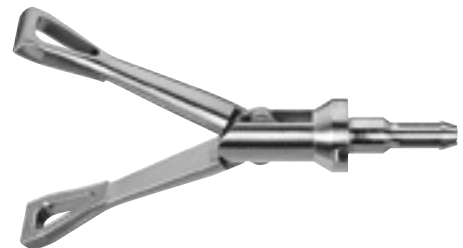

Overhold clamp / with fine teeth / curved / double action jaw

2901-263S	330 mm shaft length
2901-263	360 mm shaft length
2901-263L	450 mm shaft length



Requires shaft & handle





Integrated Needle Holder

Axial-ratchet-handle / LUER lock irrigation adapter / Ø 5 mm / 330 mm working length / four jaw types

Needle holder / with inserted TC crosswise teeth / single action jaw / straight

2901-330 straight
2901-331 angled



Needle holder / with inserted TC crosswise teeth / single action jaw / left curved

2901-334 straight
2901-335 angled



Needle holder / with inserted TC crosswise teeth / single action jaw / right curved

2901-336 straight
2901-337 angled
2901-337S angled 190 mm



Needle holder/positioner / with fixation teeth / TC dust / single action jaw / convex/concave jaw / straight

2901-340 straight
2901-341 angled





Scissors Working Part

For Ø 3 mm tube shaft instruments

Metzenbaum preparation scissors / straight / double action jaw

2901-311S 180 mm shaft length
2901-311 210 mm shaft length
2901-311L 330 mm shaft length



Metzenbaum preparation scissors / curved / double action jaw

2901-312S 180 mm shaft length
2901-312 210 mm shaft length
2901-312L 330 mm shaft length



Peritoneal scissors / straight / single action jaw

2901-313S 180 mm shaft length
2901-313 210 mm shaft length
2901-313L 330 mm shaft length



Hook scissors / single action jaw

2901-314S 180 mm shaft length
2901-314 210 mm shaft length
2901-314L 330 mm shaft length



Micro scissors / straight / single action jaw

2901-315S 180 mm shaft length
2901-315 210 mm shaft length
2901-315L 330 mm shaft length



Biopsy Cup Working Part

For Ø 3 mm tube shaft instruments

Double action jaw / without thorn

2901-111S 180 mm shaft length
2901-111 210 mm shaft length
2901-111L 330 mm shaft length



Double action jaw / with thorn

2901-113S 180 mm shaft length
2901-113 210 mm shaft length
2901-113L 330 mm shaft length





Grasping Forceps Working Part

For Ø 3 mm tube shaft instruments

Universal grasping forceps / double action jaw

2901-201S 180 mm shaft length
2901-201 210 mm shaft length
2901-201L 330 mm shaft length



Maryland grasping forceps / double action jaw

2901-202S 180 mm shaft length
2901-202 210 mm shaft length
2901-202L 330 mm shaft length



Dissecting forceps / fenestrated / double action jaw

2901-203S 180 mm shaft length
2901-203 210 mm shaft length
2901-203L 330 mm shaft length



Grasping forceps / with reservoir / blunt / double action jaw

2901-204S 180 mm shaft length
2901-204 210 mm shaft length
2901-204L 330 mm shaft length



Grasping forceps / with reservoir / sharp / double action jaw

2901-205S 180 mm shaft length
2901-205 210 mm shaft length
2901-205L 330 mm shaft length



Dolphin grasping forceps / double action jaw

2901-206S 180 mm shaft length
2901-206 210 mm shaft length
2901-206L 330 mm shaft length



Grasping forceps / sharp / crosswise teeth / 12 mm jaw length / double action jaw

2901-207S 180 mm shaft length
2901-207 210 mm shaft length
2901-207L 330 mm shaft length




Claw grasping & extracting forceps / 2x3 teeth / single action jaw

2901-208S 180 mm shaft length
2901-208 210 mm shaft length
2901-208L 330 mm shaft length




Assembling of Modular System


1

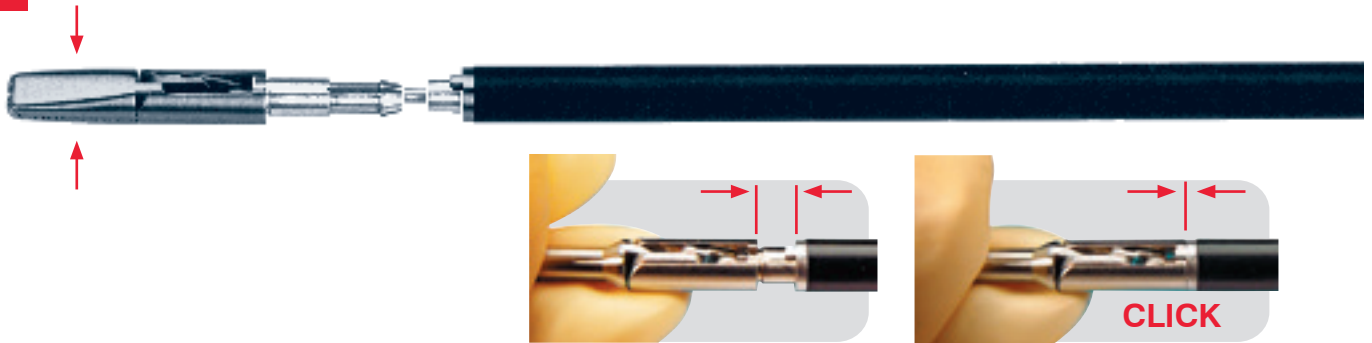

2

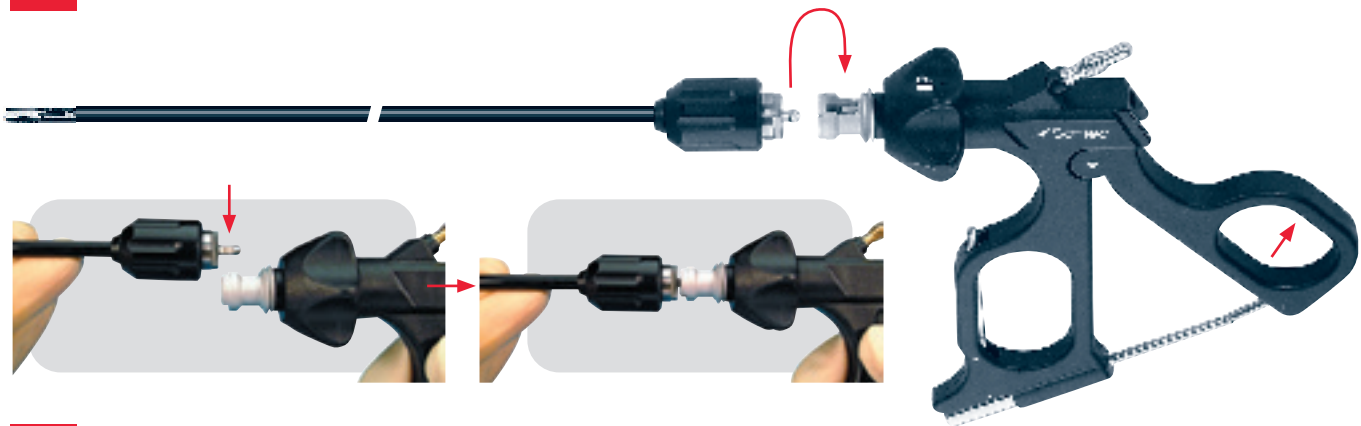

3

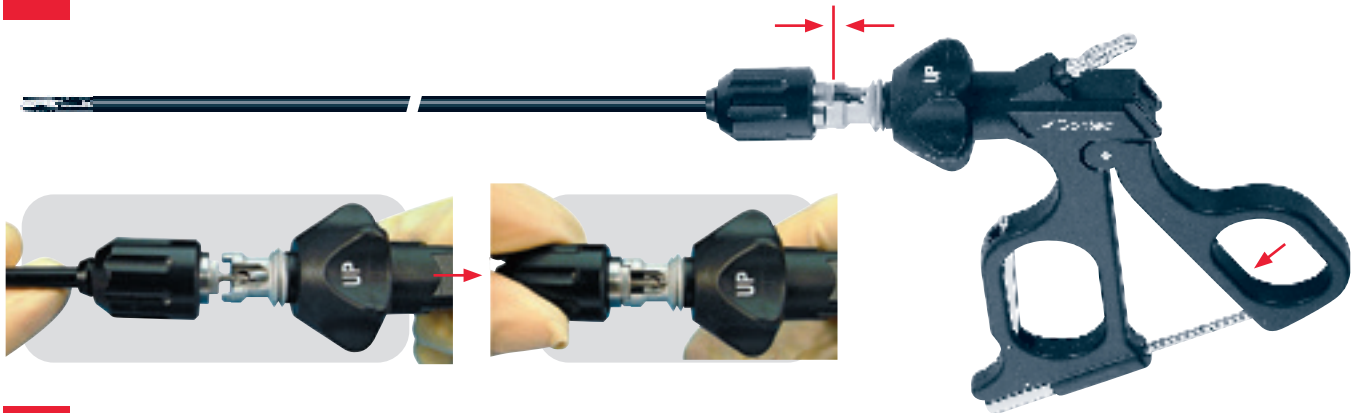

4




Disassembling of Modular System


1



2

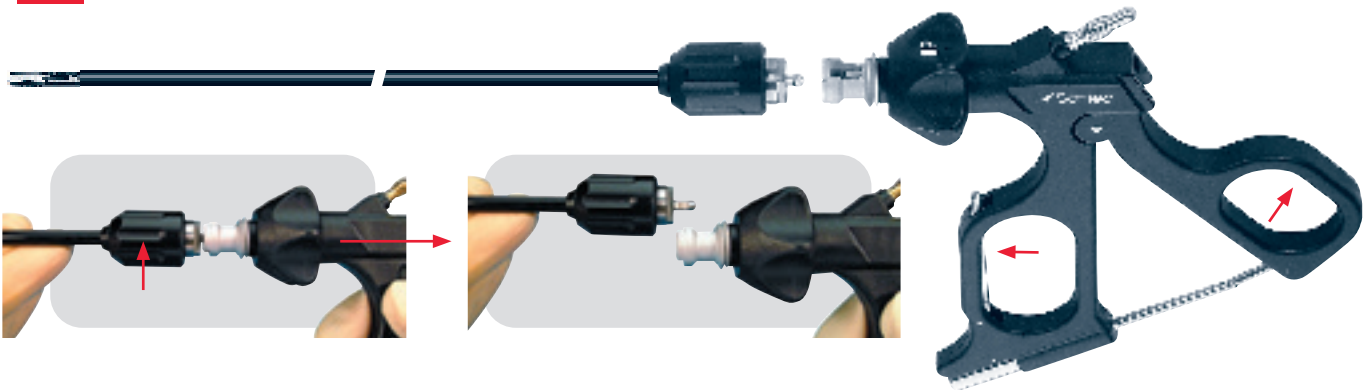


3

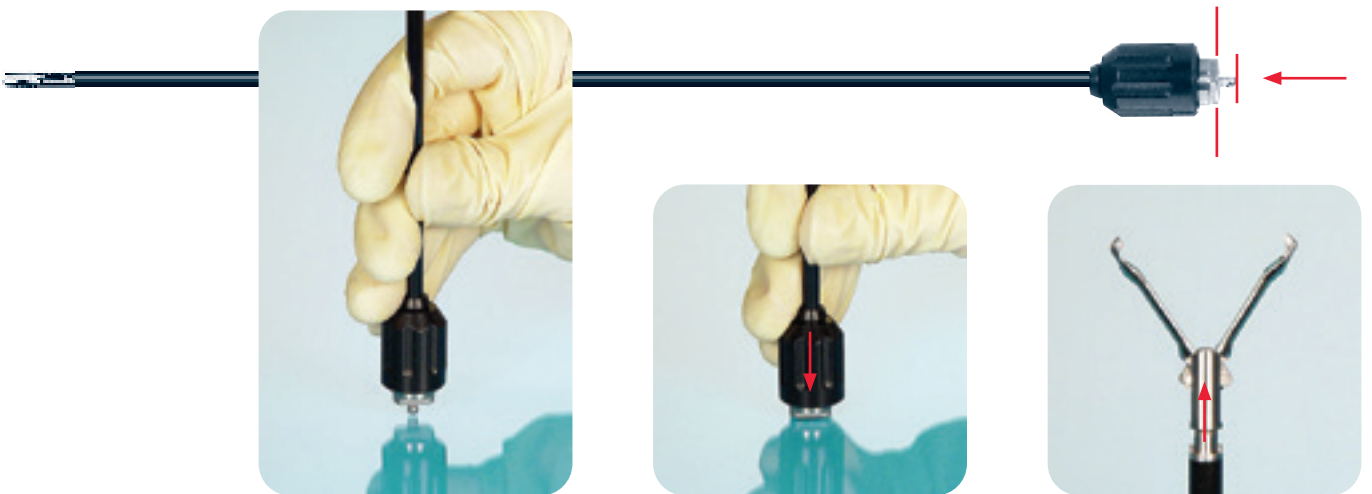

4



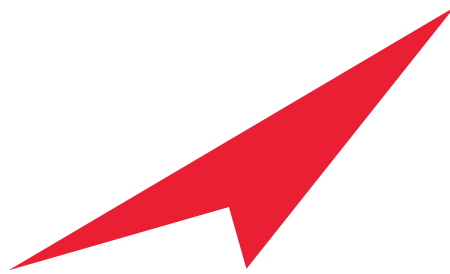

Sontec Instruments®

A TRADITION OF FINE SURGICAL INSTRUMENTS

Sontec Instruments, Inc.
7248 South Tucson Way
Centennial, Colorado 80112 USA
Toll-free / 800.821.7496
Local / 303.790.9411
Fax / 303.792.2606

SontecInstruments.com